

Instructions for installing DELTA software on your computer

This can take up to 3 days to get the licence file from Jeol.

Go to the JEOL NMR support site:


<https://nmrsupport.jeol.com/>

Scroll to the bottom of the page download the Read Me file and follow the instructions.


Firstly, you have to register as a user and after doing this you can download the processing software and install it. Please note that, if you are using a university computer, you will require an administrator password for the installation.

Once installed the software is still not licensed and so will run only for 15 minutes at a time.

Start up the Delta software and in the window that opens (Delta Console window):

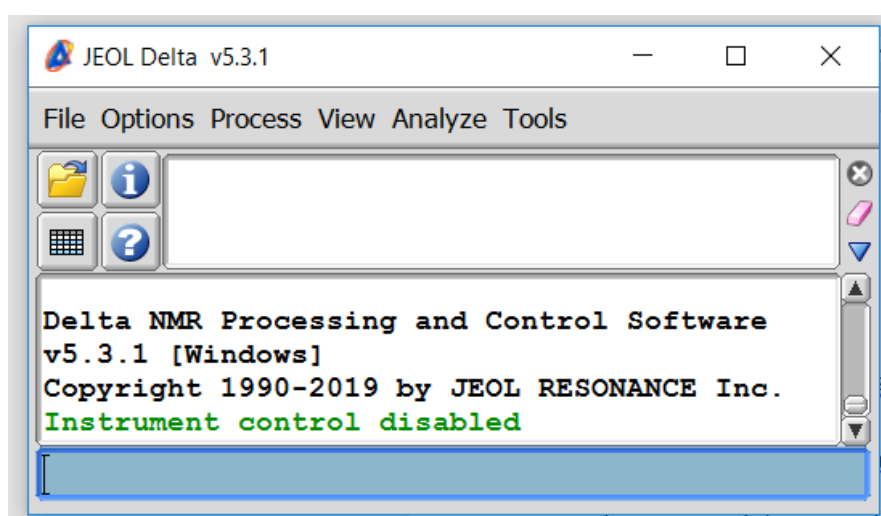

Left click on the Options tab and choose Licence Key:

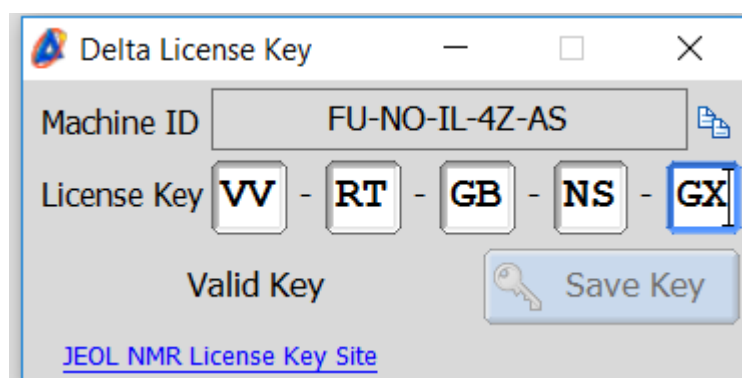

When this window opens it will show the machine ID, which is unique to your computer, but not a licence key; note the machine ID and then click on the hyperlink to the JEOL NMR License Key Site and follow the instructions, which involves entering the machine ID, to generate a license key which you then enter in the window shown here and press the Save Key softkey.

The processing software is now ready to use.

For much of this software there are multiple ways to do the same thing; we have endeavoured to suggest the easiest way(s).